

HUNTLEIGH



LymphAssist™

Garment Measuring Guide

Garment Measuring Guide


The following guide is intended as an aid to help you select the correct size garment.

The treatment type will determine the garment type for **LymphAssist™** for specialist multi chamber therapy.

When determining the correct size garment:

- Remove any bandaging, stockings or bulky dressings
- Measure the points as shown on the diagram below
 - DO NOT pull the measuring tape tight
 - Make sure the garment extends beyond the area under treatment

Refer to the garment list at the end of this document.

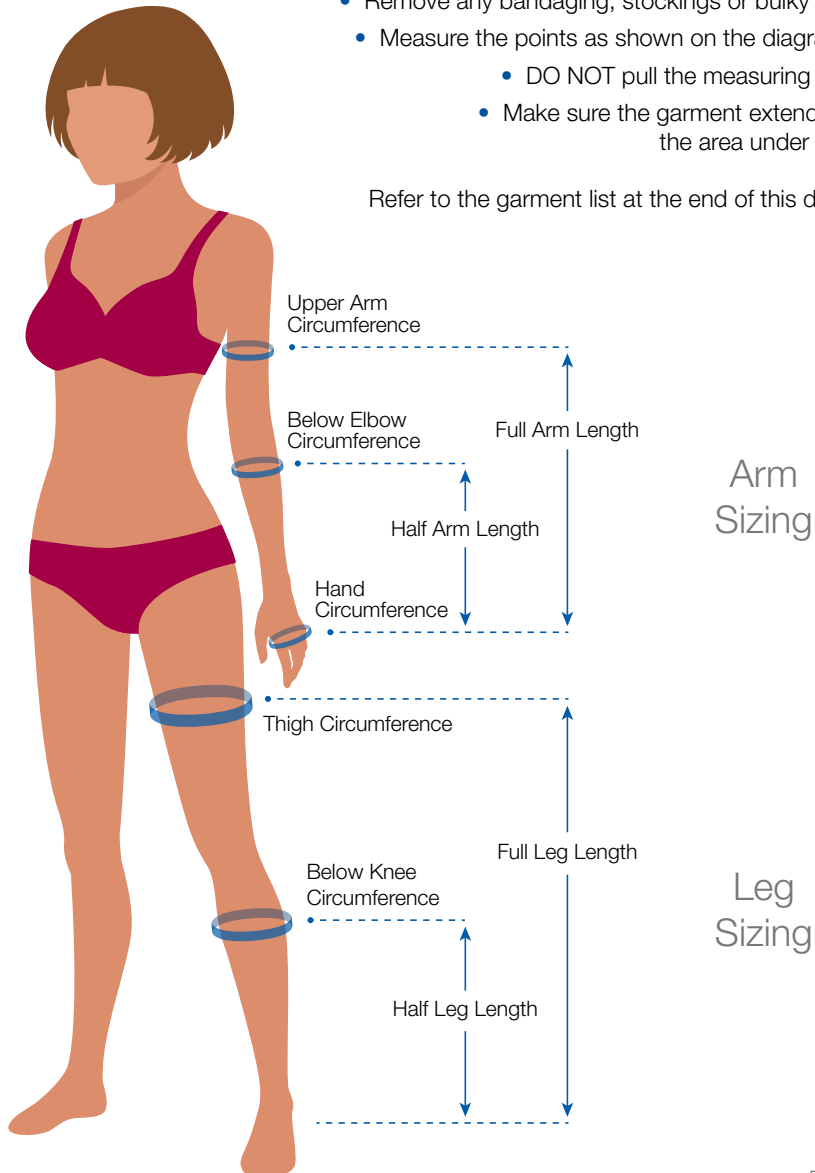

LymphAssist (Homecare/Professional) Arm

Twelve Chamber Arm Garments



Order Code	Type	Length (L)	Circumference Shoulder (H)
316A68	Full Arm	68cm	62cm
316A78	Full Arm	78cm	62cm

LymphAssist (Homecare/Professional) Arm Insert

Twelve Chamber Arm Inserts



Order Code	Type	Length (L)	Additional Width at Shoulder
316AI68	Full Arm	68cm	17cm
316AI78	Full Arm	78cm	17cm

LymphAssist (Homecare/Professional) Leg

Twelve Chamber Leg Garments



Order Code	Type	Length (L)	Thigh Circumference
316L76S	Standard Leg	76cm	71cm
316L84S	Standard Leg	84cm	71cm
316L76W	Wide Leg	76cm	79cm
316L84W	Wide Leg	84cm	79cm

LymphAssist (Homecare/Professional) Leg Insert

Twelve Chamber Leg Inserts



Order Code	Type	Length (L)	Additional Width at Thigh
316LI76	Leg Insert	76cm	19cm
316LI84	Leg Insert	84cm	19cm

HUNTLEIGH *...performance for life*

Huntleigh Healthcare Ltd. - Diagnostic Products Division

35 Portmanmoor Road, Cardiff, CF24 5HN, United Kingdom

T: +44 (0)29 20485885 **F:** +44 (0)29 20492520 **E:** sales@huntleigh-diagnostics.co.uk

W: www.huntleigh-diagnostics.com



Registered No: 942245 England & Wales. Registered Office:
ArjoHuntleigh House, Houghton Hall Business Park, Houghton Regis, Bedfordshire, LU5 5XF
©Huntleigh Healthcare Limited 2015

AN ARJOHUNTLEIGH COMPANY, MEMBER OF GETINGE GROUP

® and ™ are trademarks of Huntleigh Technology Limited

As our policy is one of continuous improvement, we reserve the right to modify designs without prior notice.

AW:1001005-1

LIT 779384/UK-1